



BEDPAN WASHERS



CDD 1050



CDD 1060



CDD 1090

Cleaning and Disinfection Devices

made by **lischka** Germany



ENGLISH VERSION

Upon reserve of technical changes

2022-02-25_Rev.00

Easy handling: To start a program, press any key on the user interface. You can also easily adjust the A0-value for disinfection (from A0-60 to A0-3000) directly on the interface. As for the program sequence the operator can follow it on the LCD display.

Hygiene & security: insulation in the chamber door reduces heat emissions and guarantees very low operating noise level, ensuring that the machine doesn't disturb patients. To protect staff from steam burn, the chamber door is only released, once the chamber has gone below a preset temperature.

Durability: The devices are made of high-grade stainless steel. This ensures more durability over the years.

Complete emptying: As per EN ISO 15883-1/-3 2009 the water tank, pumps and tubes are emptied automatically at the end of every cycle, thus avoiding any recontamination.

Cleaning and Disinfection: The practical arrangement of the multifunctional power nozzles combined with a high-capacity pump, guarantees the energy-efficient cleaning and disinfection of all contents.

PATENT

Protection: our patented recooling system protects your patients, your staff and your healthcare facility.

At the end of the disinfection phase, our patented system cools the chamber air via a condenser, so that when staff opens the chamber door, there will be no steam left inside, thus, no room pollution by aerosols.

This also means that we guide the humid warm air via a condenser as cooled air into the waste water system.

LOW OPERATING COSTS:

In addition to their durability these devices have extremely low maintenance costs. The following parts will for example not need replacement:


- No Hepa filter
- No rubber seal around chamber door

SELF DISINFECTION:

For greater hygiene and security, CDDs can self-disinfect at any preset time of the day. And if a CDD was not in use in the previous 24 hours, the device can run a complete pipe and nozzle disinfection before being used again.

Sizes, variants, all types of combinations...

Our innovative Geysir bedpan washers are available in device widths 90 cm, 60 cm and 50 cm. Depending on the model, they can be used as built-in, wall-hanging, free-standing, or under-bench. They come in all sorts of configurations (automatic doors, dosing pumps for lime descaling and/or detergent, autostart function, etc...)

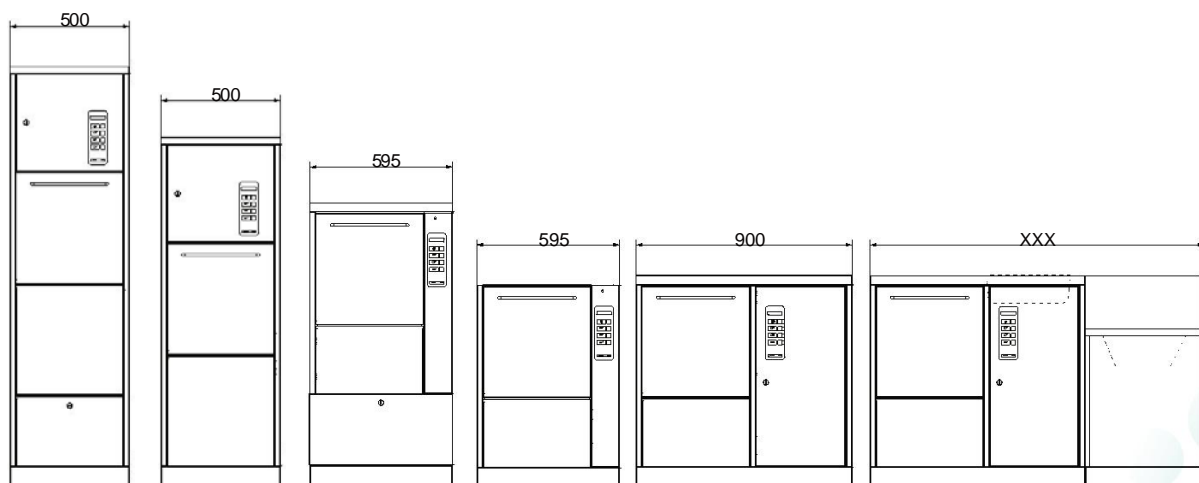

Customized solutions

We offer the perfect combination of **bedpan washers and medical furnishings**. All built in-house!

Our Geysir sluice room units combine high level hygiene with perfect practicality within sluice rooms, in hospitals, clinics or nursing homes. They are characterized by their high quality and solid construction and are available in numerous variations (material, shape, size and color).

These low maintenance, easy to care for surfaces, are ideally suited for heavy workloads and are available in finely polished stainless steel (Grain 400). or powder coated electrolytically galvanized steel sheet.

Complete Sluice rooms



Built-in CDD



Certified for...



All our devices can be equipped with special programs for **Clostridium Difficile** and Viruses (**Covid19** and **Norovirus**),

Their efficiency is tested and certified by an independent laboratory.

If you already have one of our machines, a special Kit can be installed to upgrade your device with this feature.

New
Certified programs against:

- **Clostridioides difficile**
- **Noroviruses**
- **Coronaviruses** (including Covid 19)

Available as an option on all our devices.

Certified Upgradable selectable direct on the device

Universal holding frame

Included with every device

Multifunctional

for a wide range of containers

The universal holding frame offers the highest degree of operating safety and comfort for staff using the machine.

This user-friendly frame, can accommodate the most varied utensils for everyday clinic and nursing-home life. If required, the frame can be quickly replaced, no tools needed.

Articles to be washed are emptied hygienically when the operator closes the chamber door. The operator does not have to empty the contents of the articles before placing them in the CDD. This ensures optimum hygiene, security and avoids malodors.


The ingenious flushing design of the washing chamber ensures consistent and optimal cleaning and disinfection.

Simple

to use



Easy slide-in mechanism

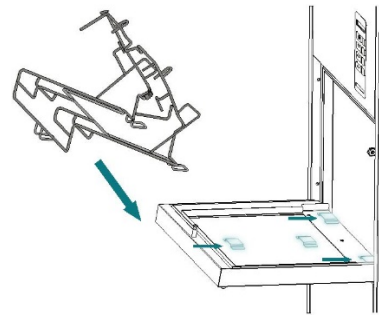

1 Bedpan and 1 or 2 Urine bottles
or 3 Urine bottles per cycle



Options

Geysir-Cut

Urine/ Secretion bag cutter

Increase the ROI of your Cleaning and Disinfection device!

With Geysir-Cut, you can disinfect urine/secretion bags and reduce your disposal waste costs.



Special Holders

Tray for small items



Basket for small parts



for suction bottles (3L)



for suction bottles (4,5L)



Our mission:

To contribute globally to the quality enhancement and affordability of infection control within the healthcare industry.



Produced with renewable energy



Lischka GmbH
Boxberger Straße 19
12681 Berlin Germany
info@made-by-lischka.de

Tel.: +49 30 981958-0
Fax.: +49 30 981958-200
www.made-by-lischka.de